

Tenutō Cool 5L316

This pharmacy refrigerator is suitable to store vaccines, medicines and laboratory products at 2°C to 8°C temperature range. Application could be found at pharmacies, pharmaceutical factories, hospitals, CDC, clinics, laboratories. Leading technologies for safer storage and higher temperature stability.

5L316
Net capacity 316L
Refrigerator 2°C to 8°C



Safer Storage Superior Temperature Uniformity and Precise Control

- » Microcomputer temperature control
- » Dual display, temperature and humidity
- » Temperature accuracy 0.1°C
- » Temperature variation $\leq 2^{\circ}\text{C}$
- » Temperature uniformity $\pm 1^{\circ}\text{C}$
- » USB port, temperature record traceable
- » Backup battery
- » Audible and visual alarm systems

Advanced Technologies Reliable and Stable Performance

- » International famous compressor
- » Top cooling system
- » DC turbo fan
- » Built-in air duct
- » LBA forming technology
- » Micro-channel condenser
- » Patented D+ heat exchange technology, high efficiency and energy saving.

Humanized Design Friendly for Daily Use

- » Electric heating glass door
- » Door self-close
- » Adjustable shelves with label holder
- » Porthole
- » Adjustable LED light
- » Wide voltage range design

Top Cooling System

Strong cooling airflow from the top to the bottom, works as an air curtain to block the external hot air from entering the refrigerator when door is opening, so as to prevent the temperature in the cabinet from rising rapidly, also contributes better temperature uniformity performance.

- » Improves temperature uniformity performance by 40%.
- » Max interior temperature difference $\pm 1^{\circ}\text{C}$
- » Prevent temperature rising rapidly when door is opening
- » Energy consumption 47% lower than similar sized normal pharmacy refrigerator.

DC Turbo Fan

- » Decrease noise level by 8%
- » Integrated DC Turbo Fan works with smaller vibration and higher efficiency
- » Safer for daily use
- » No strong electric component in chamber, safer and more reliable.

Built-In Air Duct

- » Increase the effective interior space by 10%
- » Less interior space occupancy
- » Built-in air duct plus the top cooling design, ensure temperature in the top of chamber is cool enough for storage


Dual Display

Display real-time temperature and humidity.





Micro-Channel Condenser

Compressor overheating protection, effectively improving the service life of the core part. Powerful heat dissipation for condenser, heat exchange capacity increase by more than 20% than common copper tube fin condenser of other suppliers.

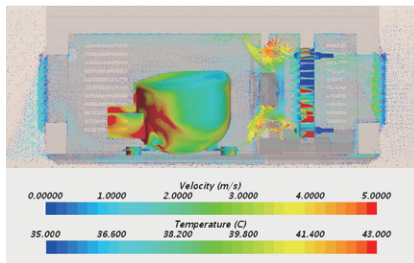

Typical Performance Characteristics at 25°C Ambient

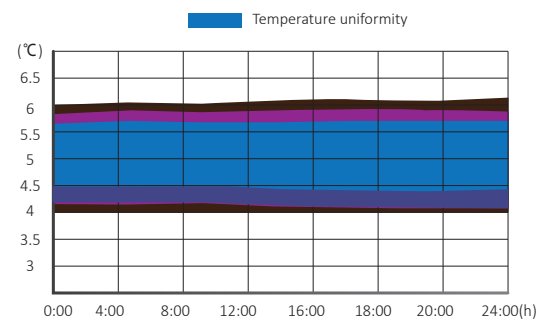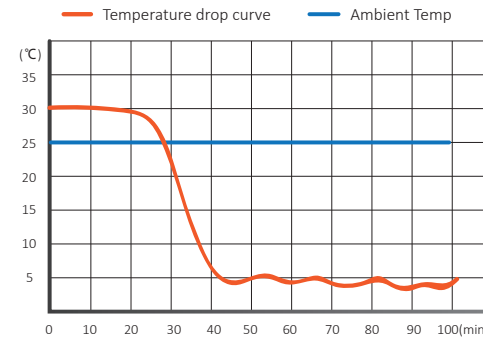

LBA Forming Technology

- » The world's leading fourth generation foaming agent.
- » More uniform density, better insulation, more energy-saving.
- » More reasonable foaming layer design, larger interior storage space.

Tenutō Cool 5L316 Specifications

Item	Model No.	5L316
Dimensions	Compartment	Refrigerator
	Load Capacity (40HQ)	71
	Fridge Storage Volume (L)	316
	Product Dimensions WxDxH (mm)	595*603*1920
	Interior Dimensions WxDxH (mm)	500*464*1343
	Net/Gross Weight	78/90
Technical Data	Compartment	Refrigerator
	Cooling Type	Forced air cooling
	Defrost Type	No frost
	Climate Class	N
	Refrigerant/Amount (g)	R600/36g
	Compressor Type	Fixed frequency
	Sensor Quantity	6
Product Performance	Compartment	Refrigerator
	Temperature Range	2°C to 8°C
	Temperature Uniformity/°C	2
	Temperature Volatility/°C	2
	Temperature Accuracy	0.1
Control	Compartment	Refrigerator
	Controller	Microprocessor
	Display	LED (Temperature and Humidity)

Item	Model No.	5L316
Electrical Data	Compartment	Refrigerator
	Rated Current(A)	1.5
	Rated Power(W)	190
Alarms	Compartment	Refrigerator
	Power Failure Alarm	Y
	Sensor Failure Alarm	Y
	Environment Temperature Alarm	Y
	Door Ajar Alarm	Y
	High Temperature Alarm	Y
	Low Temperature Alarm	Y
	Remote Alarm	Optional
Accessories	Compartment	Refrigerator
	Shelves/Drawers Quantity	5
	Anti-Condensation Door	Electrical heating
	Self-Lock	Y
	Label Card Holder	Y
	Cable Port	Y
	Castors/Type	4/Omni-directional
	USB	Y
	RS-485 Port	Y
	Printer	Optional
	Backup Battery	Y
	Backup Battery Endurance	24 hour





**For more information, visit:
swisslog-healthcare.com/tenuto**

